



Lutera

The 5-Minute Home Vision Check

Your bonus gift — find your baseline today, then re-check it every 30 days.

What's inside

- **Test 1** — Eye Strain Score: how hard modern life is working your eyes
- **Test 2** — Amsler Grid: a classic check of your central vision
- **Test 3** — Near Vision Chart: the smallest print you can read
- **Test 4** — Contrast Check: the faintest letters you can see
- Your 90-Day Tracker: log today's results and re-test as you go

Before you start

- Print this at 100% ("Actual size") — never "Fit to page."
- Check the ruler below: the gold bar must measure exactly 10 cm.
- Sit in bright, even light. Wear the glasses you normally read with.
- Test one eye at a time — cover the other with your palm (don't press).
- Keep this sheet. You'll re-test on days 30, 60, and 90.



PRINT CHECK — this bar should measure exactly 10 cm (100 mm) with a ruler

Test 1 — Eye Strain Score

Check YES for anything you've noticed in the last month. Takes 60 seconds.

- ☐ 1. My eyes feel tired or heavy by the evening.
- ☐ 2. My eyes feel dry, gritty, or burning after screen time.
- ☐ 3. Things look blurry for a moment when I look up from my phone.
- ☐ 4. I zoom in or enlarge text on my phone to read it.
- ☐ 5. I squint to read menus or labels in dim light.
- ☐ 6. I get headaches after reading or computer work.
- ☐ 7. Headlights at night seem to glare or halo more than they used to.
- ☐ 8. I need brighter light to read comfortably than I used to.
- ☐ 9. I catch myself rubbing my eyes during the day.
- ☐ 10. I hold my phone or books farther away to focus.
- ☐ 11. My eyes take longer to adjust when I walk into a dark room.
- ☐ 12. I spend more than 6 hours a day looking at screens.

Your score = number of YES boxes: _____ / 12

0–3 Low strain — nice. Protect it.


4–7 Moderate strain — your eyes are asking for backup.

8–12 High strain — make the next three tests your baseline and re-check monthly.

Test 2 — The Amsler Grid

A classic central-vision check used in eye clinics for over 75 years.

- Hold this page 30 cm (12 inches) from your face. Wear your reading glasses.
- Cover your LEFT eye. Stare only at the gold dot in the center.
- Without moving your eye off the dot, notice the lines around it.
- All lines should look straight, evenly spaced, with no blurry or missing patches.
- Repeat with your RIGHT eye covered. Note anything odd in your tracker.

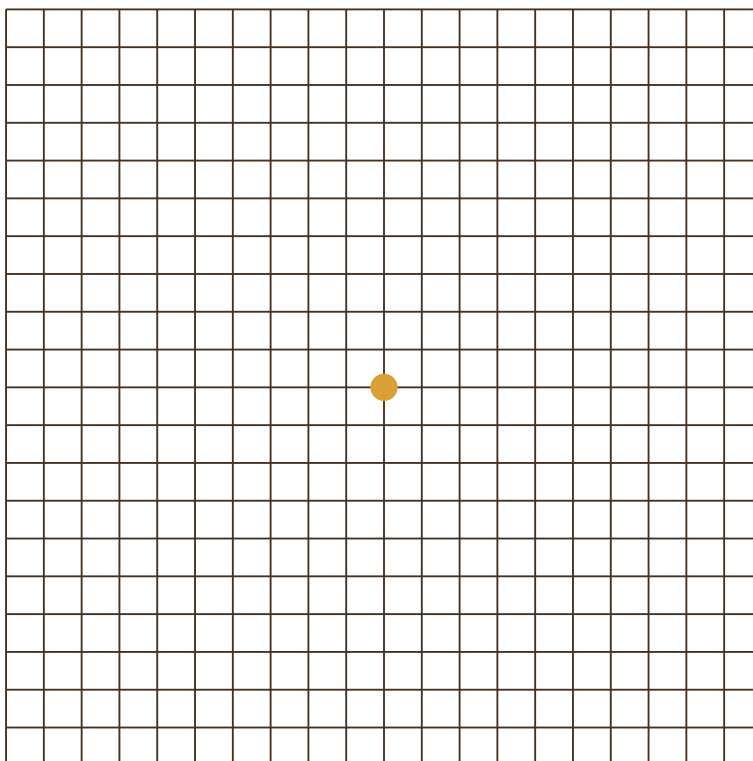

IMPORTANT: if any lines look wavy, broken, blurry, or missing — with either eye — book an appointment with an eye-care professional promptly. Don't wait for day 30.

Test 3 — Near Vision Chart

Hold at 40 cm (16 in), one eye at a time, reading glasses ON. Letter sizes are true to scale.

20/200	LOOK UP	large headline
20/100	REST YOUR EYES	sub-headline
20/70	BLINK MORE OFTEN	children's book print
20/50	YOUR EYES WORK HARD	newspaper print
20/40	EVERY SINGLE WAKING HOUR	magazine print
20/30	GIVE THEM WHAT THEY NEED	paperback print
20/25	PROTECT THE VIEW AHEAD	medication label
20/20	SEE WELL FOR LIFE	fine print

Your result = the smallest line you can read comfortably, per eye.

Write it in your tracker (e.g. "R: 20/30, L: 20/40"). Struggling two lines apart between eyes, or well short of 20/40 with your usual glasses? That's worth an eye exam.

Test 4 — Contrast Check

Real life is low-contrast: dusk, fog, headlights, grey roads. Hold at 40 cm, one eye at a time.

- Each row fades. Read across, row by row, until you can no longer make out the letters.
- Your result = the last row you read correctly. Note it per eye in your tracker.
- Printer inks vary — always compare against your OWN earlier print-outs, not anyone else's.

ROW 1	D	N	R	Z	K
ROW 2	H	C	V	O	S
ROW 3	K	Z	N	D	R
ROW 4	O	S	H	C	V
ROW 5	R	D	K	Z	N
ROW 6	V	H	S	O	C

Fading contrast sensitivity often shows up YEARS before a regular eye chart catches anything — it's one of the first things people notice as night driving gets harder.

Your 90-Day Tracker

Your eyes didn't get this way overnight — and change is gradual. Measure, don't guess.

Fill in the Baseline row today. Re-run all four tests on days 30, 60, and 90 — same room, same light, same glasses, and this same printed sheet.

	Date	Strain score (/12)	Amsler notes (R / L)	Smallest line (R / L)	Faintest row (R / L)
Baseline					
Day 30					
Day 60					
Day 90					

While you track

- Take Lutera daily — consistency beats intensity. Same time every morning is easiest.
- Follow 20-20-20: every 20 minutes, look at something 20 feet away for 20 seconds.
- Blink fully and often at screens. Most people's blink rate halves while scrolling.
- Get outside in daylight — and wear sunglasses when it's bright.

The fine print (please actually read this bit)

This booklet is a self-screening and tracking tool for general wellness purposes only. It is not a medical device, not an eye examination, and cannot diagnose, treat, or rule out any condition. It does not replace a comprehensive exam by an optometrist or ophthalmologist — most adults should have one regularly, whatever their results here. If you notice sudden vision changes, distortion on the Amsler grid, flashes, new floaters, or eye pain, seek professional care immediately.

Lutera is a dietary supplement. These statements have not been evaluated by the Food and Drug Administration. This product is not intended to diagnose, treat, cure, or prevent any disease.